


minerva™
The Uterine Health Company

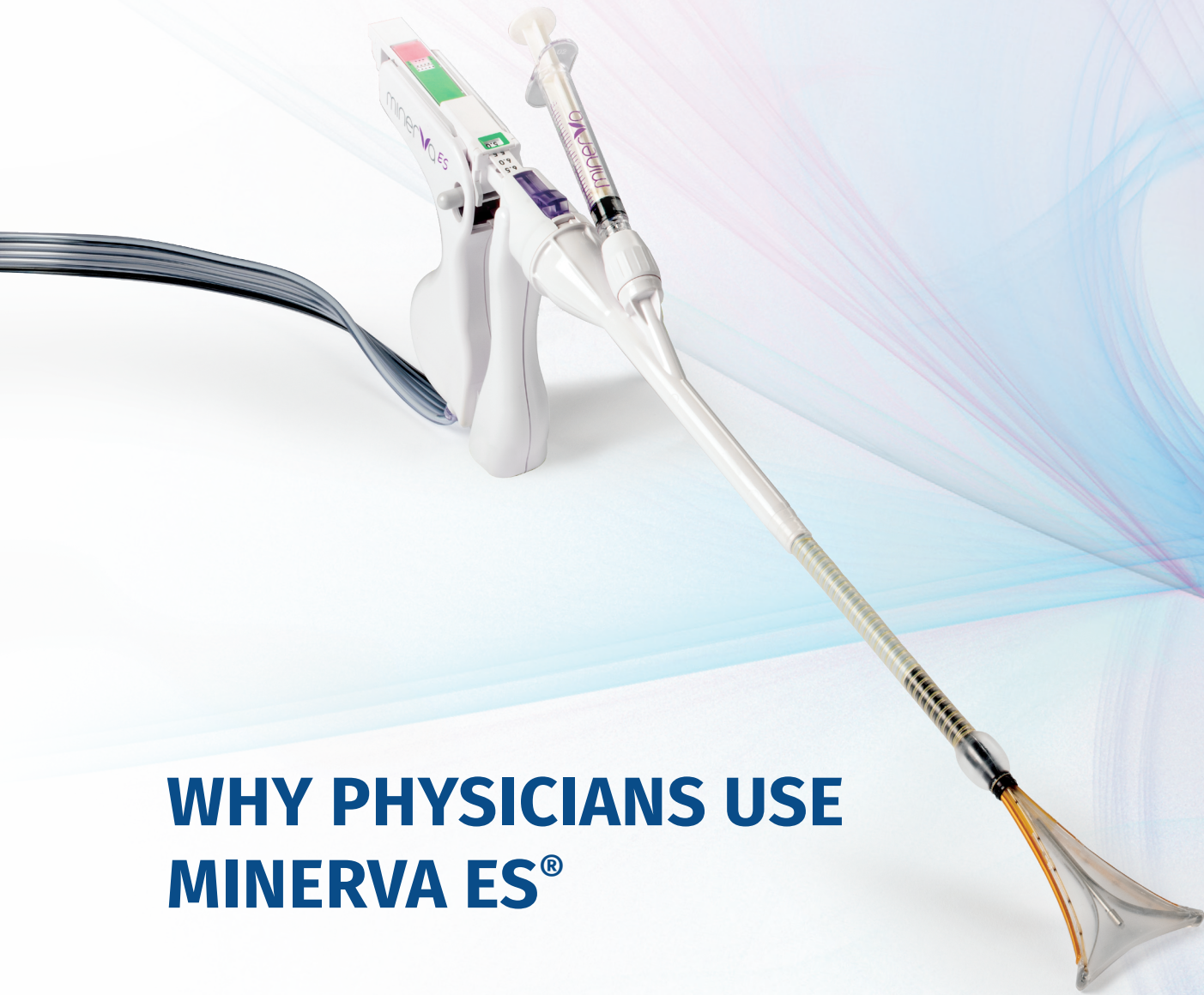

WHY PHYSICIANS USE MINERVA ES®

Minerva ES®
Endometrial Ablation System

Mastery made easy


Soft, silicone array allows easy insertion, deployment and removal.

- ✓ No coarse mesh - array does not stick to tissue

- ✓ Eliminates difficult seating


Minerva ES Cervical Seal Balloon produces a tight cervical seal, reducing procedural delays.

- ✓ Hands-free cervical seal
- ✓ Designed to seal the cervix at the internal os and protect from unintended thermal effects

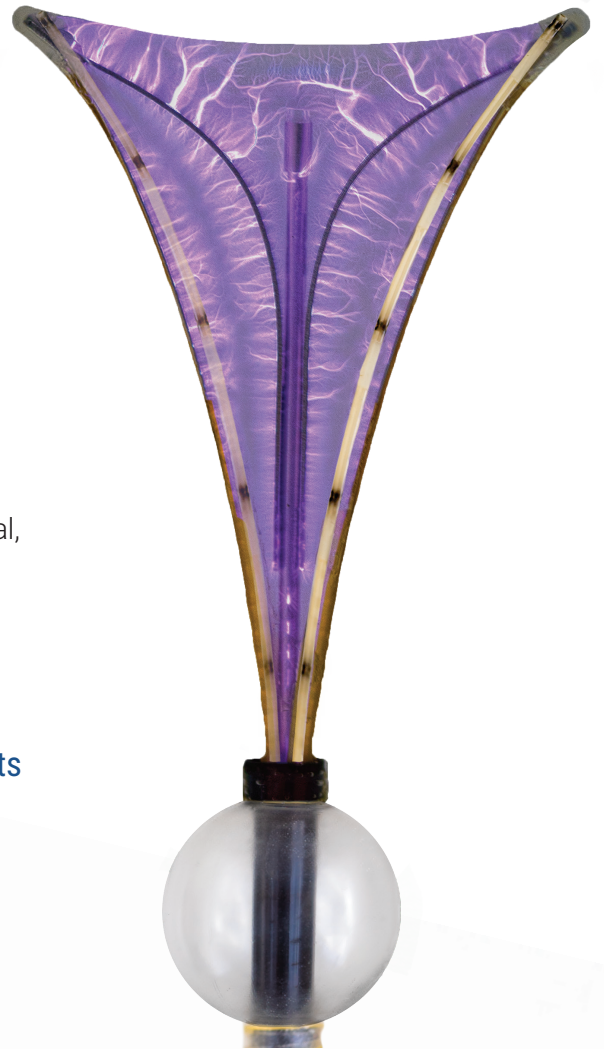

Cervical Seal Indicator confirms a complete seal prior to Uterine Integrity Test (UIT).

- ✓ Indicates absence of CO2 leak
- ✓ Streamlines the treatment and saves time

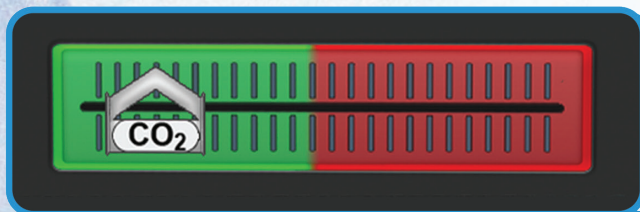


Safety enhanced

Patented CO2 Extension Tubes extend along the entire length of the array.

- ✓ CO2 flows from lower uterine segment to the tips of the array via 29 gas delivery ports
- ✓ Multi-port CO2 delivery design allows for even distribution of gas, increasing the chance of perforation detection



Patented two-stage Uterine Integrity Test improves detection of perforations.

- ✓ Each stage passes when CO2 flow remains below the threshold for 10 seconds
- ✓ Competitive one-stage system must maintain pressure threshold for only 4 seconds

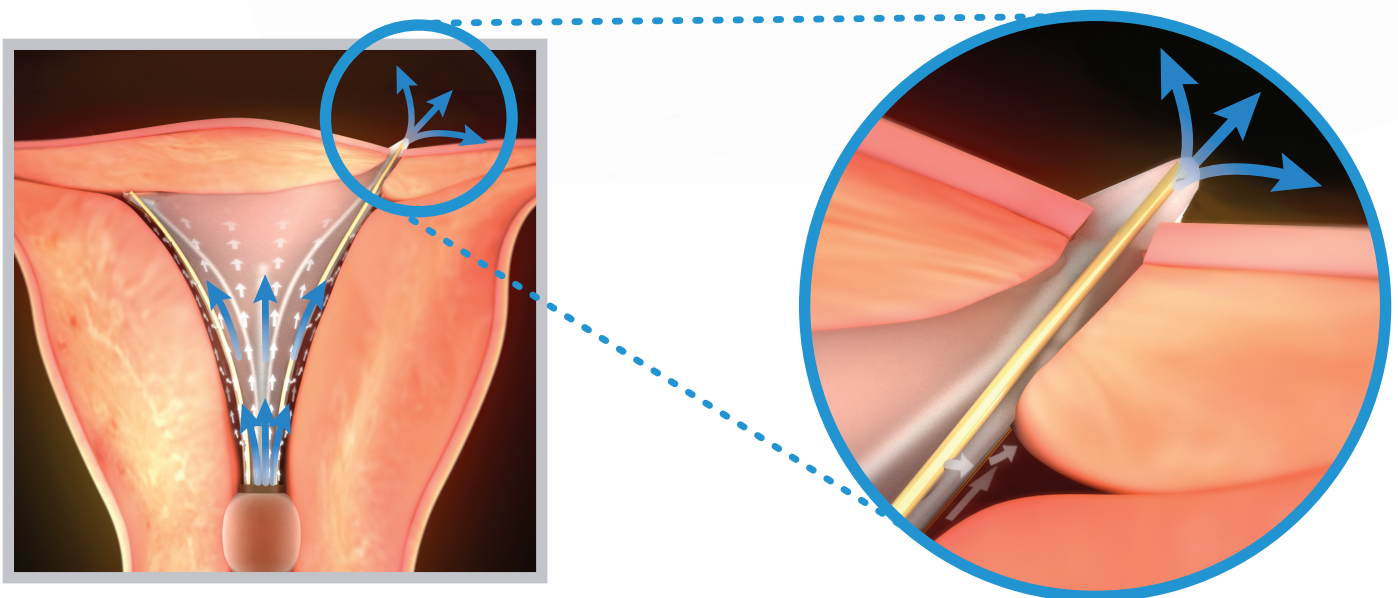

Clinical data

Significant reduction in post-ablation hysterectomy rate

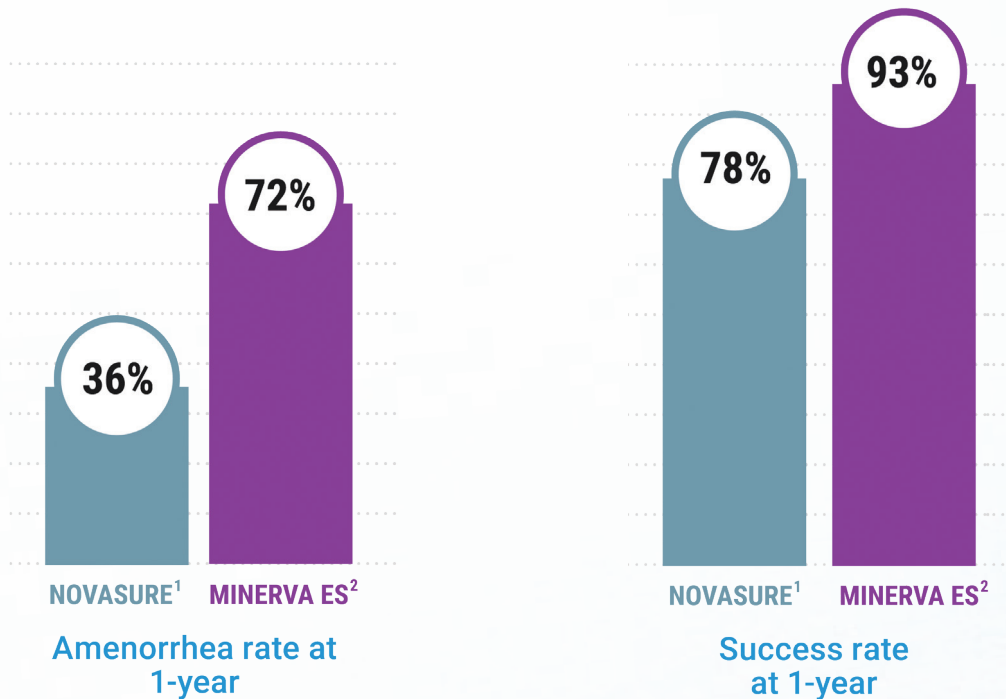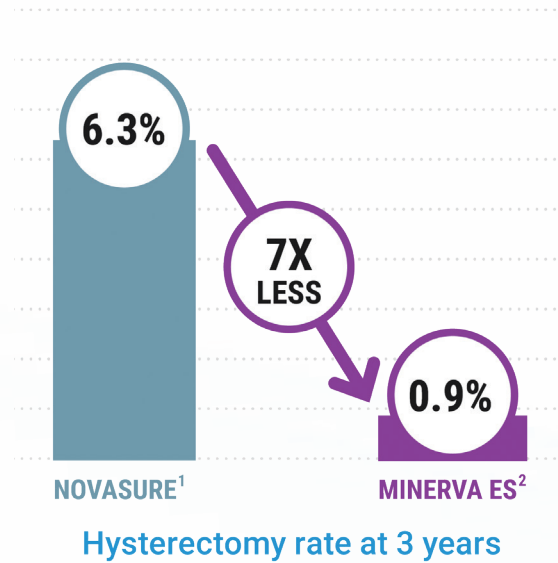

1. FDA Approved Labeling: NovaSure Impedance Controlled Endometrial Ablation System [Instructions For Use and Controller Operator's Manual]. Marlborough, MA: Hologic, Inc.; 2014.

2. FDA Approved Labeling: Minerva Endometrial Ablation System [Operator's Manual]. Santa Clara, CA: Minerva Surgical Inc.; L0120, 2020.

Order Minerva ES and experience the ease and efficacy of a modern ablation device

Contact your Minerva sales representative or place your order directly:

Call: 855-646-7874 Fax: 866-465-2875

Email: customerservice@minervasurgical.com

| Description | Part Number |
|------------------------------------|-------------|
| Minerva ES Handpiece (3-pack) | MIN3PAK |
| Controller | MIN180S |
| CO ₂ Canisters (5-pack) | MINC02C |
| Argon Canisters (5-pack) | MINARGC |

Additional resources

Visit the [Resource Center](#) at minervasurgical.com for more information on our complete suite of devices for the treatment of AUB.



The most advanced FDA approved endometrial ablation system with best-in-class clinical results. PlasmaSense™ technology delivered with dynamic power is designed to ensure complete coverage and consistent depth of ablation.

Learn more at www.minervasurgical.com

Minerva Surgical
4255 Burton Dr.
Santa Clara, CA 95054

© 2020 Minerva Surgical
or its affiliates. All rights reserved.

K0099 Rev. B

minerva™
The Uterine Health Company

Federal law restricts this device to sale by or on order of a physician.

Minerva ES® is a registered trademark of Minerva Surgical, Inc. All other trademarks are the property of their respective owners.